

Within this section, we will provide some common usage examples of SWATH-Lib.

## 0.1 Creating a library containing Y-ion for N-glycosylated peptide

It has been known that *Y*-ion based library can be used for analysis of glycopeptide without prior glycan knowledge.[**RN474**] Here we have a precursor sequence with the id **test\_sequence** and the aa sequence **TTENDTFWKEF**, or in fasta format

```
>test_sequence
TTENDTFWKEF
```

Upon input into the fasta content box and process, display similar to Figure 0.1 should appear.

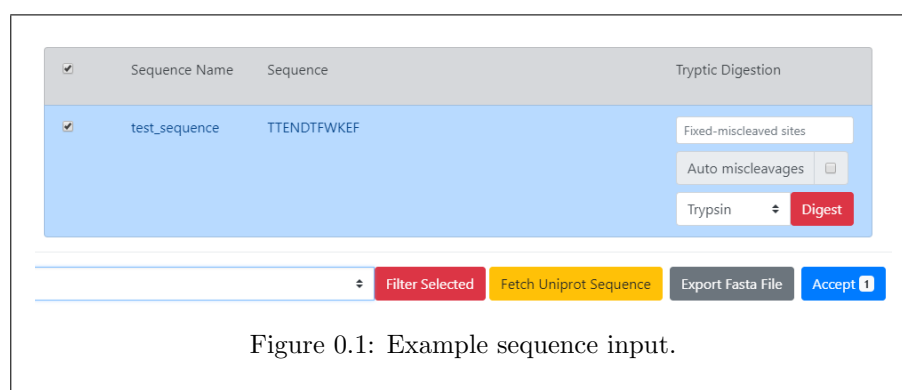


Figure 0.1: Example sequence input.

This sequence has an N-glycosylation sequon at position 4 and a tryptic digestion site right after position 9. With this library, we aimed to generate all HexNAc Y0, Y1, and Y2 transitions for the sequence at RT of 10 across all SWATH windows. The max charge allowed would be 2 and the annotated sequence at the modification position would include both RT and SWATH window value.

After in-silico digestion with trypsin, we would obtain two sequences (Figure 0.2), one before and one after the digestion site. The name of the sequences have also been automatically modified to show the starting and stopping positions of each digested fragment. Using the filter rule for N-glycosylation sequon, the GUI would identify and select all sequences with a sequon within those already selected.

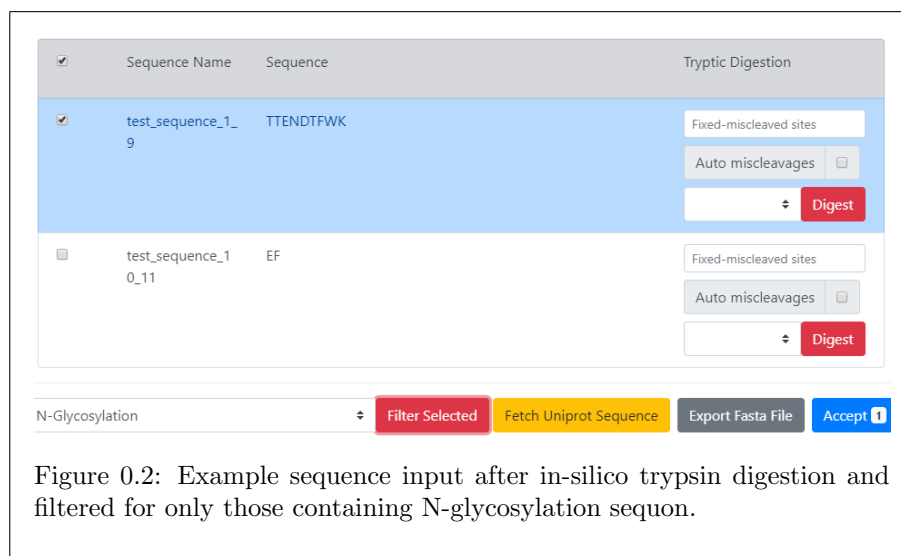


Figure 0.2: Example sequence input after in-silico trypsin digestion and filtered for only those containing N-glycosylation sequon.

We currently have a tryptic digested sequence containing N-glycosylation sequon but without any modification and experiment settings applied to it. These settings can be selected within the **Queryset Input Settings** panel.

In Figure 0.3, we have the following settings selected,

- HexNAc Y0, Y1, and Y2 modifications.
- Retention time at 10.
- All SWATH windows.
- Y fragmentation ion-type
- Output sequence format include RT and SWATH window values at variable modification sites.

**Queryset Input Settings**

**Static Modifications**

- Propionamide (71.0371 Da)
- Carbamidomethyl (57.02 Da)

**Non-Ytype Variable Modifications**

- Phosphorylation (79.9663 Da)
- Sulfation (79.9568 Da)
- Carboxylation (43.9898 Da)
- Methylated Carboxylation (58.0055 Da)

**Ytype Variable Modifications**

- O-Mannose (Y1) (162.05 Da)
- HexNAc (Y0) (0 Da)
- HexNAc (Y1) (203.0794 Da)
- HexNAc (Y2) (406.1587 Da)

**Retention Time**

- 10
- 11
- 12
- 13

**SWATH Windows**

- 1149 - 1175
- 1174 - 1200
- 1199 - 1225
- 1224 - 1250

**Extra Mass**

0

**Max Charge**

2

**Precursor Charge**

2

**Oxonium Ions**

- Hex-2H<sub>2</sub>O (127.0389 m/z)
- HexNAc-(H<sub>2</sub>O)<sub>2</sub>-COH<sub>2</sub> (138.055 m/z)
- Hex-H<sub>2</sub>O (145.0495 m/z)
- Hex (163.0601 m/z)

**Fragmentation ion-type**

Y

**Output Sequence Format At Variable Modifications**

Retention Time & SWATH

**Apply Modifications**

Figure 0.3: Example modification settings for N-glycosylation library with 3 HexNAc Ytype transitions at RT 10 and across all default SWATH windows.

After application of modification settings, the result would be similar to Figure 0.4. Multiple HexNAc transitions are placed at N4 and the residue is colored green.

**test\_sequence\_1\_9** Edit

1 2 3 4 5 6 7 8 9 10 11 12 13 14 15 16 17 18 19 20

1- T T E N D T F W K

☐ b,y-run non-Ytype mods only

☐ Generate only for oxonium ion

**HexNAc Y0: 1**

☐ Filled all positions

☐ Multiple Pattern

**HexNAc Y1: 1**

☐ Filled all positions

☐ Multiple Pattern

**HexNAc Y2: 1**

☐ Filled all positions

☐ Multiple Pattern

Y

Figure 0.4: Example of a sequence with the settings above applied.

On submission, the progress bar would start. When a query is finished, the bar under it would turn green. The moment all queries result have finished and been collected, a **Save Compiled Results** button would appear under the **Submit Queries** button. The user can click on the button to save the library in .txt tabulated file format.

## Creating a byY-library

For a library with both *by* and *Y* transitions, the user might want to avoid overlapping of the transition *m/z*, especially *y* and *Y* transitions. The **b,y-run non-Ytype mods only** option can be enabled for the query. This option would disable Ytype modification mass within calculation of transition *m/z* for *by*-series.

## 0.2 Adding oxonium ions to Y-library

In MS data of glycosylated peptides, Oxonium ions are fragmentation of the modification block during ionization process. The presence of oxonium ions fragment within the spectral graph of a peptide fragment could provide additional evidence for the peptide to be glycosylation and potential knowledge on the glycan composition. With SWATHLib, using the query settings interface, you can add these oxonium ions to your library of interests to look for potential glycosylation evidence.

Let start again with the example sequence from the section above.

```
>test_sequence
TTENDTFWKEF
```

Following similar input procedure as of Figure 0.3 we would obtain the digested peptide containing at least one N-glycosylation sequon. However, different to what depicted previously, our settings would now include the oxonium ions Hex, HexNAc, and HexHexNAc.

The screenshot displays the SWATHLib query settings interface with the following configurations:

- Static Modifications:** Propionamide (71.0371 Da), Carbamidomethyl (57.02 Da).
- Non-Ytype Variable Modifications:** Phosphorylation (79.9663 Da), Sulfation (79.9568 Da), Carboxylation (43.9898 Da), Methylated Carboxylation (58.0055 Da).
- Ytype Variable Modifications:** O-Mannose (Y1) (162.05 Da), HexNAc (Y0) (0 Da), HexNAc (Y1) (203.0794 Da), HexNAc (Y2) (406.1587 Da).
- Retention Time:** 10, 11, 12, 13.
- SWATH Windows:** 1149 - 1175, 1174 - 1200, 1199 - 1225, 1224 - 1250.
- Extra Mass:** 0.
- Max Charge:** 2.
- Precursor Charge:** 2.
- Oxonium Ions:** HexNAc-H2O (186.0761 m/z), **HexNAc (204.0866 m/z)**, NeuAc-H2O (274.0921 m/z), NeuAc (292.1027 m/z).
- Fragmentation ion-type:** (Empty field).
- Output Sequence Format At Variable Modifications:** SWATH Windows.

Figure 0.5: Example of a sequence with the oxonium ions selected. By selecting these settings, every single modified variation of the peptide at each selected window and retention time would also be accompanied by entries of selected oxonium ions.

These settings together with the query can now be submitted for process to produce a library with Y-ion together with those selected oxonium ions.

## Oxoniome library

An oxoniome library is one that containing only oxonium ions. It can be used to quickly gain an overview of glycosylation within complex samples. SWATHLib can be used to create an oxoniome library from peptide of interests.

The screenshot shows the SWATHLib interface for a peptide sequence named 'test\_sequence\_1\_9'. The sequence is displayed as '1 - T T E N D T F W K' with positions 1 through 20 indicated above. The 'N' at position 4 is highlighted in green. To the right of the sequence, there are three modification entries, each with a green label: 'HexNAc Y0:1', 'HexNAc Y1:1', and 'HexNAc Y2:1'. Each entry has a 'Filled all positions' checkbox and a 'Multiple Pattern' checkbox. The 'Generate only for oxonium ion' checkbox is checked. Below these, there is a 'Selected oxonium ions:' section with a list of options: 'Hex', 'HexNAc', and 'HexHexNAc'. At the bottom right, there is an 'Ion Type' input field and a 'Clear All Modifications' button.

Figure 0.6: Example of a sequence with the oxonium ions selected and is set for only oxonium generation.

In Figure 0.6, the option **Generate only for oxonium ion** was checked to instruct the back-end only generating oxonium ions entry for each modification pattern without calculating any *Y* or *by*-transitions. Currently this setting is only available on an individual sequence basic. The sequence and all modifications input would only be used to identify different modification pattern of the input peptides.

SWATHLib can also generate oxoniome library for arbitrary peptides as long as no modifications and transitions are put on the arbitrary peptides in the settings.